

Specification

Chambers	1-3
Mode	4
Time Range	1 min - 99 hours or continuous
Pressure Range	30 - 150mmHg
Hold Time	2s (Default)
Interval Time	10s (Default)
Interface	2.2" OLED screen
Rechargeable Battery	Yes (Lithium battery: 3500mAh)
Automatic Memory	Yes
Dimension	270(L) x 148(W) x 145(H)mm
Weight	2.2kg

Certification



CE(MDD)



ISO13485 / MDSAP



Health Canada

Longest
FURNISH & BEHOLD SERVICES

Guangzhou Longest Medical Technology Co., Ltd.

Address: 301&401 of Building 2& Building 3, No.96, Chuangqiang Road, Ningxi Street, Zengcheng District, Guangzhou, Guangdong Province, 511399, China

Tel: 86-20-66353999 Fax: 86-20-66353920

E-mail: export@longest.cn Website: www.longestmedical.com



LGT-2201DVT

Prophylaxis of Deep Vein Thrombosis (DVT)

Health Life Longest Care

Longest

LGT-2201DVT

LGT-2201DVT is a special pressotherapy system with gradient pressure that is designed to decrease the risk of DVT (Deep Venous Thrombosis) and Pulmonary Embolism (PE) before, during and after surgery.

1 Professional Overlapping DVT Sleeve

⊗ Disposable Sleeve

(Foot+Calf / Thigh+Calf / Foot / Calf / Thigh)

⊕ Reusable Sleeve




2 4 Compression Modes


MODE 1




MODE 2

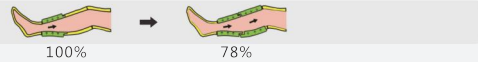

MODE 3

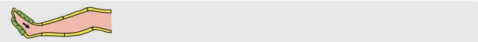

MODE 4

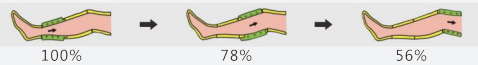

3 Gradient Pressure

Ensures the patient achieves optimal therapeutic and comfort during treatment.

4 Real-time Pressure Monitoring

The real-time pressure monitoring system can automatically fill the pressure when the pressure is lower than the setting and stop after reaching the setting.

5 Bedside / Backside Hook



6 Flexible Treatment Time

1min-99h or continuous output

7 DC / AC Supply

8 Automatic Memory for Last Setting

9 Screen Lock Function

If you do not operate the device for about 10s, it will automatically lock to prevent misoperation.

10 Quick Connector

Securely and easily snaps into the hose

